

KESYON? RELE (800) 232-1700

# SIDEKICK

KÒMANSE AK  
SMARTWATCH ALÈT MEDIKAL OU A



# Chaje aparèy ou an

## Pou pi bon rezilta chaje Sidekick ou an chak swa.



Yon fwa ou resevwa Sidekick ou a, ploge aparèy ki pou chaje li a nan yon priz. Mete Sidekick la nan aparèy ki pou chaje li a pandan bouton èd SOS la bay fas li sou devan. Asire w ke Sidekick ou a chita byen fèm nan aparèy ki pou chaje li a pou li ka kòmanse chaje. Lè lap chaje kòrèkteman **Limyè WOUJ** ki sou aparèy ki pou chaje li a ap limen epi Sidekick la ap anonse "Charging."

Lè Sidekick la konplètman chaje, limyè **WOUJ** ki sou aparèy ki pou chaje li a ap etèn.



Asire w ke fil ki nan dèyè aparèy ki pou chaje li a byen ranje.

KESYON?

RELE (800) 232-1700

# Konnen aparèy ou an

# Bouton



Ou kapab benyen ak **Sidekick la** men fòk ou pa plonje li nan dlo.



Peze **Bouton Silans** lan de fwa pou dezaktive/aktive alèt vokal.

# Limyè

# Tès Aparèy ou a



## LIMYÈ STATI

Limyè VÈ **a limen etèn** lè Sidekick la ap fonksyone nòmalman.

## LIMYÈ LOKALIZASYON

Limyè BLE **a limen etèn** lè sèvis lokalizasyon yo disponib.

## LIMYÈ POWER

Ou dwe mete Sidekick ou a chaje lè limyè **WOUJ** la ap limen etèn.



## BOUTON TÈS

Pou teste Sidekick ou a, peze epi kenbe bouton tè s la ki arebò Sidekick ou a. Ou pral konekte ak sant tè s nou an.



## REMAK

Pou pi bon rezilta, teste **Sidekick** ou a chak mwa.

# Nan Ka Ijans

**Si w bezwen asistans, peze epi kenbe bouton èd la jiskaske ou santi ou yon vibration oswa ou tande tonalite aktivasyon yo.**



## **ENPÒTAN**

**Pou anile yon aktivasyon ki fèt pa aksidan, peze epi lage bouton èd la lè yo mandew fè sa.**



## **Deteksyon Tonbe** (opsyonèl)

Pou aparèy ki gen Deteksyon Tonbe ki aktive sou li, si ou pa kapab peze bouton èd la epi li detekte ou sot tonbe, aparèy la ap voye yon alèt ijans pou tonbe a bay yon Espesyalis Swen.

Pou pi bon rezilta, mete bouton èd ou a nan yon ti kòd epi pa deyò rad ou.



Deteksyon tonbe a pa detekte 100% tonbe yo. Si yo kapab, itilizatè yo ta dwe toujou peze bouton an lè yo bezwen èd.

# Etèn/Limen

Sidekick ou a dwe toujou rete limen.  
Si ou bezwen **ETÈN** Sidekick ou a peze  
epi **KENBE** bouton Tès ak èd yo ansanm  
jiskaske ou santi yon vibrasyon.

Tout limyè ki sou aparèy la ta dwe etèn  
nan yon minit.

Pou **RELIMEN Sidekick la** ou ka peze bouton  
Tès ak bouton èd la ansanm, oswa mete  
aparèy la nan lòt aparèy ki pou chaje li a.

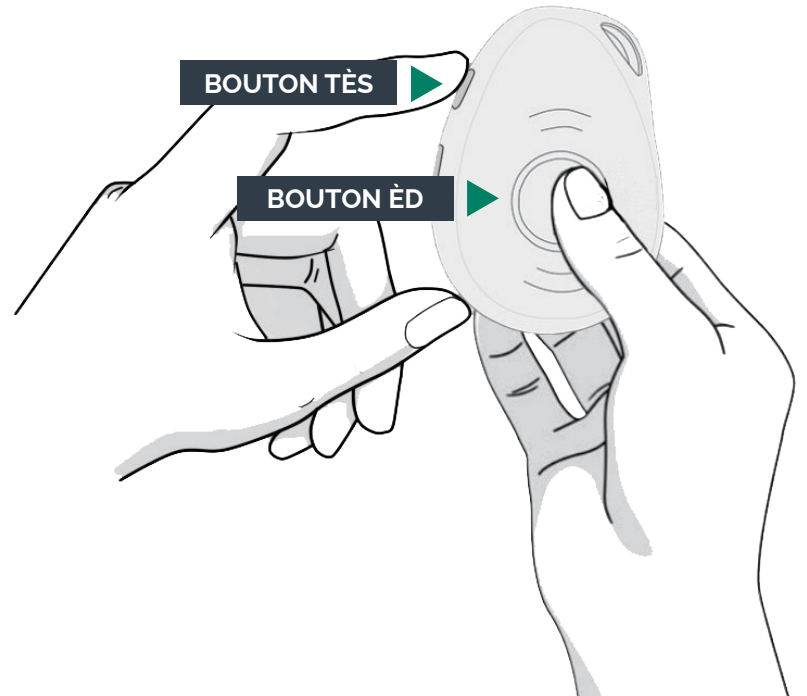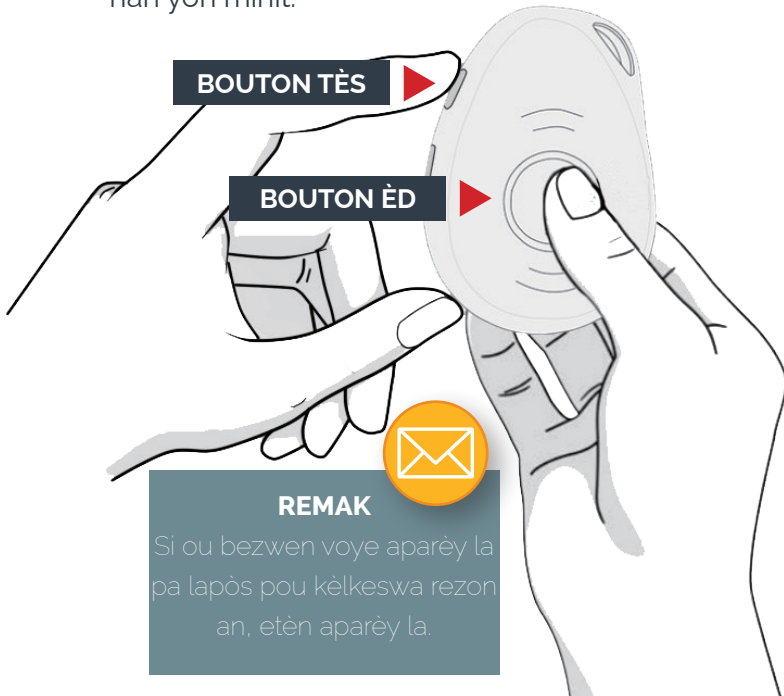

KESYON?

RELE (800) 232-1700

**Nou la pou nou ede.  
Rele (800) 232-1700  
pou asistans.**

Aparèy sa a depann de kouvèti selilè pou li fonksyone.  
Nan kote kouvèti selilè a limite, ou ka remake kapasite aparèy la diminye.